


RA-b-23

## Automated\_Total RNA Extraction from Liver of Mouse



Depending on sample and storage conditions, nucleic acid may not be extractable. Therefore, we cannot guarantee accurate data. The extracted nucleic acid contains unintended acid (ex: when extracting DNA, RNA is also extracted). KURABO

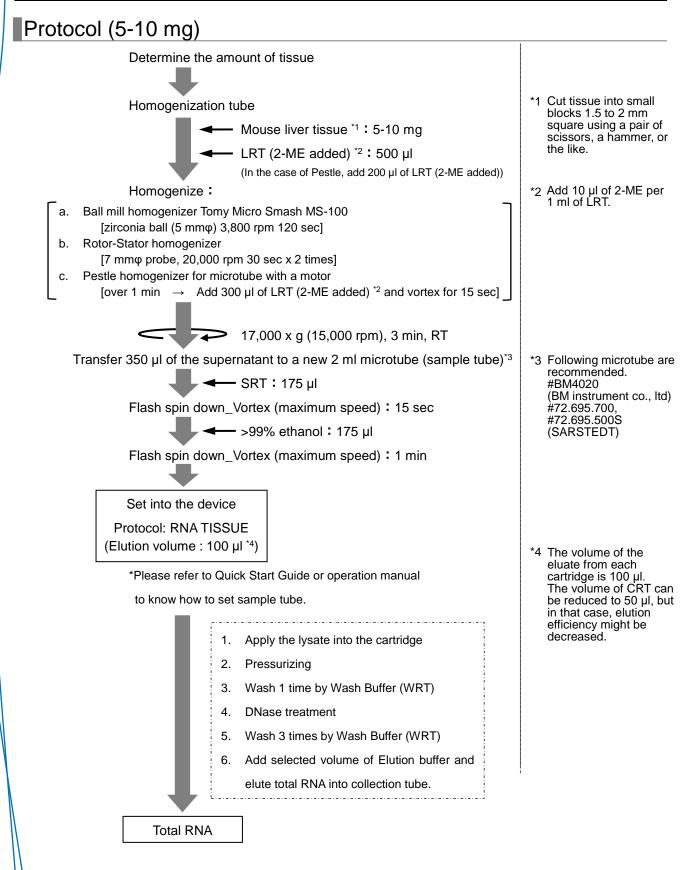

Results

## The yield of total RNA

Sample ID	#1	#2	#3	#4	Average	S.D.
Yield (µg)	35.0	29.8	29.9	28.1	30.7	3.01

## Protein contamination: A260/280 /Chaotropic salt contamination: A260/230

Sample ID	#1	#2	#3	#4	Average
A260/280	2.36	2.32	2.33	2.30	2.33
A260/230	2.15	1.90	2.16	2.05	2.07

## Common protocol is usable for the following

Mouse testis, Mouse Liver, Mouse Brain, Mouse Lung, Mouse Spleen

